## Precision

### **Water Baths - Shaking, Circulating, Coliform and General Purpose**

Manual Part Number U01316, November 5, 2015

#### **Multilingual Quick Start Guides**

**Installation**

**Operation**

**Preventive Maintenance**

**Troubleshooting**

# Table of Contents

## Quick Start Guides

|  |            |
|--|------------|
| <b>Preface</b>                         | <b>P-1</b> |
| Compliance .....                       | P-1        |
| Unpacking .....                        | P-1        |
| After-Sale Support .....               | P-1        |
| Feedback.....                          | P-2        |
| <b>Chapter 1 Safety</b>                | <b>1-1</b> |
| Safety Factors.....                    | 1-1        |
| <b>Chapter 2 General Information</b>   | <b>2-1</b> |
| Description .....                      | 2-1        |
| Intended Use .....                     | 2-1        |
| Specifications .....                   | 2-2        |
| <b>Chapter 3 Installation</b>          | <b>3-1</b> |
| Bath Installation .....                | 3-1        |
| Ventilation.....                       | 3-1        |
| Cover Installation .....               | 3-1        |
| Electrical Requirements.....           | 3-2        |
| Approved Fluids .....                  | 3-3        |
| Filling Requirements.....              | 3-6        |
| Filling with Thermal Beads.....        | 3-6        |
| Draining .....                         | 3-7        |
| Accessories .....                      | 3-8        |
| Drain Kit Accessory .....              | 3-8        |
| <b>Chapter 4 Operation</b>             | <b>4-1</b> |
| Controller.....                        | 4-1        |
| Start Up.....                          | 4-2        |
| Viewing/Changing the Shaker Speed..... | 4-3        |
| Changing the Set Point.....            | 4-3        |
| Selecting a Preset Set Point .....     | 4-3        |
| Changing a Setting .....               | 4-4        |
| Selecting a Preset.....                | 4-6        |
| Changing a Preset .....                | 4-6        |
| Temperature Display .....              | 4-7        |
| Shut Down.....                         | 4-7        |

|                  |                                    |            |
|------------------|------------------------------------|------------|
| <b>Chapter 5</b> | <b>Preventive Maintenance.....</b> | <b>5-1</b> |
|                  | Pump Inlet Filter.....             | 5-1        |
|                  | Cleaning.....                      | 5-1        |
|                  | Electrical Power Cord.....         | 5-2        |
|                  | Optical Disc.....                  | 5-2        |
|                  | Shaker Assembly.....               | 5-2        |
| <b>Chapter 6</b> | <b>Troubleshooting.....</b>        | <b>6-1</b> |
|                  | Error Displays.....                | 6-1        |
|                  | Checklist.....                     | 6-2        |
| <b>Warranty</b>  |                                    |            |

# Precision Water Bath Quick Start Guide

This quick start guide is intended for initial start up only. For all other procedures you must refer to the manual. Also, if any of these steps are not clear refer to the manual before proceeding.

## Filling (Only approved fluids are Filtered/Single Distilled Water and De-ionized Water)

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To avoid spilling, place your samples/trays into the bath before filling.

Slowly fill the reservoir. When adding, point the opening of a container away from yourself.

For Shaking Water Baths, fill the reservoir with a minimum of one inch of water and the maximum level lower than the tray shaft fitting. Make allowance for water splashing from tray oscillations.

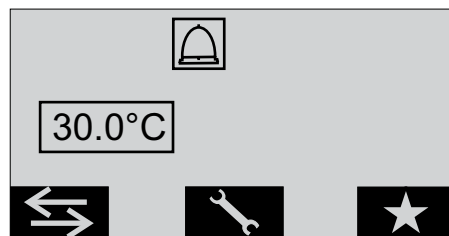
For Circulating and Coliform Baths, fill the reservoir with a minimum of two inches of water, ensuring the suction screen is covered, and a maximum level one inch lower than the upper tank surface.

For GP Baths, fill the reservoir with a minimum of one inch of water and a maximum level one inch lower than the tank upper surface.

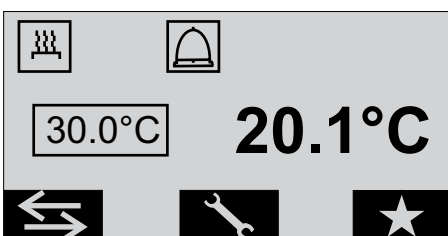
## Starting

Do not run the bath until fluid is added to the reservoir. Have extra fluid on hand. If the bath does not start refer to the manual.

- Place the circuit protector located on the rear of the bath to the I position.
- The screen will momentarily display **Thermo** SCIENTIFIC and then the home screen appears displaying the set point and the audible alarm status. The set point is the desired reservoir fluid temperature.

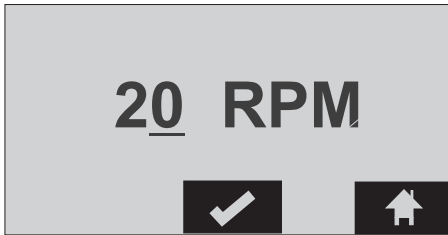
 appears on Shaking Water Baths only.

Press  to start the bath and display the fluid temperature. If the reservoir fluid temperature is below the set point the heater  will start.

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You can view/change the speed value with the bath running or not. From the home screen press the soft key below  to display:

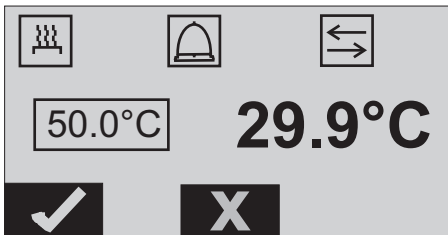
Press, and hold, the arrow keys to bring up the desired value. The values are 0 and 30 to 200. Once the desired value is displayed press  or  to save the change and return to the home screen.

## Changing the Set Point (Shaker, Circulating and GP Baths only)

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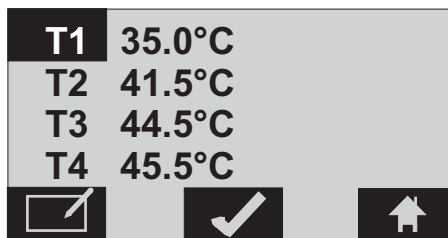
You can change the set point with the bath running or not. From the home screen press either arrow key to display:

Press, and hold, the arrow keys to bring up the desired set point value. The range is 5°C to 100°C. Once the desired value is displayed, in this case 50.0°C, press the soft key below  to save the change or press the soft key below  to not save the change. The display returns to the home screen.

**Note** After a 30 second delay if neither key is pressed the display returns to the home screen, any change is not saved.

## Selecting a Preset

You can select a preset with the bath running or not.

From the home screen press  :

Use the arrow keys to select the desired preset. Once the desired preset is highlighted press the soft key below  to save the change and return to the home screen.

**Note** After a 30 second delay if neither key is pressed the display returns to the home screen, any change is not saved.

From the home screen you can cycle through the set points by pressing the arrow keys (Coliform baths only).

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# Preface

## Compliance

The Declaration of Conformity is available upon request.

## Unpacking

The bath is supplied with an electrical power cord. Do not discard the packaging until the cord is located and the bath is operating.

If the bath shows external or internal damage contact the transportation company and file a damage claim. Under ICC regulations, this is your responsibility.

**The bath does not have handles. Take into account its weight when unpacking and transporting. We recommend two people lift heavier baths from the bottom.**

## After-sale Support

Thermo Fisher Scientific is committed to customer service both during and after the sale. If you have questions concerning the bath operation, or questions concerning spare parts or Service Contracts, call our Sales, Service and Customer Support phone number, see next page for contact information.

Before calling, please obtain the following information:

- bath model number
- bath serial number
- power source voltage

The model and serial number are located on its nameplate label on the rear of the bath.

## Feedback

We appreciate any feedback you can give us on this manual. Please e-mail us at:

[tcmanuals@thermofisher.com](mailto:tcmanuals@thermofisher.com)

Please include the manual part number and the revision date listed on the front cover.

### **Thermo Fisher Scientific**

25 Nimble Hill Road

Newington, NH 03801

Tel : (800) 258-0830 or

(603) 436-9444

Fax : (603) 436-8411

[www.thermoscientific.com/tc](http://www.thermoscientific.com/tc)

### **Sales, Service, and Customer Support**

25 Nimble Hill Road

Newington, NH 03801

Tel: (800) 258-0830

Sales: 8:00 am to 5:00 pm

Service and Support: 8:00 am to 6:00 pm

Monday through Friday (Eastern Time)

Fax: (603) 436-8411

[service.tc.us@thermofisher.com](mailto:service.tc.us@thermofisher.com)

# Chapter 1 Safety

## Safety Factors

Make sure you read and understand all instructions and safety precautions listed in this manual before installing or operating your bath. If you have any questions concerning the operation of your bath or the information in this manual, please contact us.

There are no special personal protective equipment requirements needed to perform normal operation. We do recommend always wearing eye protection and gloves.

**DANGER** indicates an imminently hazardous situation which, if not avoided, will result in death or serious injury.

**WARNING** indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

**CAUTION** indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. It is also used to alert against unsafe practices.

The lightning flash with arrow symbol, within an equilateral triangle, is intended to alert the user to the presence of non-insulated "dangerous voltage" within the bath's enclosure. The voltage magnitude is significant enough to constitute a risk of electrical shock.

This label indicates the presence of hot surfaces.

This label indicates read the manual.

This label indicates a hand crushing hazard. (Shaking Water Bath only)

**Observe all warning labels.**

**Never remove warning labels.**

**The bath's construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection will not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided.**

**Do not use the bath as a sterile or patient connected device. In addition, the bath is not designed for use in Class I, II or III locations as defined by the National Electrical Code.**

**The circuit protector located on the rear of the bath is not intended to act as a disconnecting means. The bath's power cord is used as the disconnecting device, it must be easily accessible at all times.**

**Never operate the bath with a damaged line cord.**

**Never place the bath in a location or atmosphere where excessive heat, moisture, or corrosive materials are present.**

**Never operate the bath without fluid in the reservoir. The only approved fluids are filtered/single distilled water and de-ionized water. General Purpose Baths may also use thermal beads.**

**Never operate the bath or add fluid to the reservoir with panels removed.**

**Operating the baths at high temperature will cause condensation on the underside of the cover. Remove the cover with care so the water falls back into the reservoir and not on the bath surface.**

**When operating at high temperatures do not touch the cover surface, always use the cover's knob/handle.**

**The user is responsible for any hazard rising from heated material.**

**Do not clean the bath with solvents, only use a soft cloth and water.**

**Drain the bath before it is moved. Drain the bath before it is transported and/or stored. Store the bath in the temperature range -25°C to 60°C (with packaging), and <80% relative humidity.**

**Always turn the bath off and disconnect the supply voltage from its power source before moving the bath or before performing any service or maintenance procedures.**

**Transport the bath with care. Sudden jolts or drops can damage its components.**

**An operating shaking bath has moving components. Do not insert appendages between the reservoir wall and rack. We recommend stopping the shaker table when loading/removing samples or vessels.**

**Refer service and repairs to a qualified technician.**

**Performance of installation, operation, or maintenance procedures other than the ones described in this manual may result in a hazardous situation and will void the manufacturer's warranty.**

# 2

## Chapter 2 General Information

### Description

All Thermo Scientific™ Precision Water Baths provide constant temperature to applications placed in the reservoir. All have digital displays, four programmable set point temperatures, acoustic and optical alarms, and offer adjustable high temperature protection.

The interior of the bath is constructed of stainless steel and is designed for operation with water. The body is made from galvanized steel and is painted for added protection.

The nameplate label on the bath identifies its electrical requirements.

### Intended Use

The baths are intended for use in research and quality control.

They are intended for use by a qualified lab technician trained in basic laboratory procedures and safety protocols. Also, for indoor use in a laboratory environment on a bench top only.

## Shaking Water Bath Specifications

|  | SWB 15/SWB 15S/DUB 15                   | SWB 27                                  |
|--|---|---|
| <b>Reservoir Fluid Control Temperature</b> °C<br>°F          | Ambient +5 to 100<br>Ambient +9 to 212  |   |
| <b>Set Point Temperature Range</b> °C<br>°F                  | +5 to 100<br>+9 to 212                  |   |
| <b>Ambient Temperature Range</b> °C<br>°F                    | +15 to 45<br>+59 to 113                 |   |
| <b>Stability/Uniformity</b><br>@ 37°C                        | ±0.1/±0.05                              |   |
| <b>Heating Output</b> watts                                  | 1200                                    | 1500                                    |
| <b>Bath Volume</b> liters                                    | 15                                      | 27                                      |
| <b>Overall Bath Dimensions</b><br>(L x W x H) cm<br>inches   | 39.4 x 63.2 x 24.9<br>15.5 x 24.9 x 9.8 | 39.4 x 93.8 x 24.9<br>15.5 x 36.9 x 9.8 |
| <b>Bath Work Area Dimensions</b><br>(L x W x H) cm<br>inches | 29.2 x 30.5 x 16.5<br>11.5 x 12.0 x 6.5 | 29.2 x 61.0 x 16.5<br>11.5 x 24.0 x 6.5 |
| <b>Tray Dimensions</b><br>(L x W x H) cm<br>inches           | 28.9 x 30.8 x 14.8<br>11.3 x 12.1 x 5.8 | 28.9 x 61.5 x 14.8<br>11.3 x 24.2 x 5.8 |
| <b>Tray Depth</b> cm<br>inches                               | 8.9<br>3.5                              |   |
| <b>Shaking Speed</b> rpm                                     | 0 and 30 - 200                          |   |
| <b>Stroke Length</b> mm                                      | 26                                      |   |
| <b>Max Tray Load</b> kg<br>lb                                | 9<br>20                                 | 14<br>30                                |
| <b>Approximate Weight</b> kg<br>lb                           | 22<br>48                                | 27<br>62                                |
| <b>Electrical Requirements</b> (VAC/Hz)<br>(Voltage ±10%)    | 100-115/50-60<br>or<br>200-230/50-60    |   |
| <b>Compliance</b>  | CE RoHS<br>UL WEEE                      |   |
| <b>Maximum Relative Humidity</b><br>(Non Condensing)         | 80% (up to 31°C)<br>80% (up to 88°F)    |   |
| <b>Operating Altitude</b> meters<br>feet                     | Sea Level to 2000<br>Sea Level to 6560  |   |
| <b>Overvoltage Category</b>                                  | II                                      |   |
| <b>Pollution Degree</b>                                      | 2                                       |   |
| <b>Storage Temperature Range</b> °C<br>°F                    | -25 to +60<br>-13 to +140               |   |

Low-end temperatures require supplemental cooling.

Specifications obtained at sea level using water.

Bath depth includes bezel, height does not include lid.

Thermo Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Thermo Fisher Scientific reserves the right to change specifications without notice.

## Circulating Bath Specifications

|  | CIR 19                                  | CIR 35                                  | CIR 89                                    |
|--|---|---|---|
| <b>Reservoir Fluid Control Temperature</b> °C<br>°F          | Ambient +5 to 100<br>Ambient +9 to 212  |   |   |
| <b>Set Point Temperature Range</b> °C<br>°F                  | +5 to 100<br>+9 to 212                  |   |   |
| <b>Ambient Temperature Range</b> °C<br>°F                    | +15 to 45<br>+59 to 113                 |   |   |
| <b>Stability @ 37°C</b>                                      | ±0.1                                    |   |   |
| <b>Uniformity @ 37°C</b>                                     | ±0.05                                   | ±0.05                                   | ±0.1                                      |
| <b>Heating Output</b> watts                                  | 1200                                    | 1500                                    | 1500                                      |
| <b>Bath Volume</b> liters                                    | 19                                      | 35                                      | 89  |
| <b>Overall Bath Dimensions</b><br>(L x W x H) cm<br>inches   | 39.4 x 63.2 x 24.9<br>15.5 x 24.9 x 9.8 | 39.4 x 93.8 x 24.9<br>15.5 x 36.9 x 9.8 | 54.6 x 116.0 x 30.0<br>21.5 x 45.7 x 11.8 |
| <b>Bath Work Area Dimensions</b><br>(L x W x H) cm<br>inches | 30.5 x 38.7 x 19.2<br>12.0 x 15.3 x 7.6 | 30.5 x 69.2 x 19.2<br>12.0 x 27.3 x 7.6 | 48.3 x 91.4 x 24.1<br>19.0 x 36.0 x 9.5   |
| <b>Approximate Weight</b> kg<br>lb                           | 20<br>45                                | 31<br>68                                | 38<br>85                                  |
| <b>Electrical Requirements</b> (VAC/Hz)<br>(Voltage ±10%)    | 100-115/50-60<br>or<br>200-230/50-60    |   |   |
| <b>Pump Pressure</b> mbar                                    | 262                                     |   |   |
| <b>Pump Flow</b> lpm   | 10.6                                    |   |   |
| <b>Compliance</b>  | CE      RoHS<br>UL      WEEE            |   |   |
| <b>Maximum Relative Humidity</b><br>(Non Condensing)         | 80% (up to 31°C)<br>80% (up to 88°F)    |   |   |
| <b>Operating Altitude</b> meters<br>feet                     | Sea Level to 2000<br>Sea Level to 6560  |   |   |
| <b>Overvoltage Category</b>                                  | II                                      |   |   |
| <b>Pollution Degree</b>                                      | 2                                       |   |   |
| <b>Storage Temperature Range</b> °C<br>°F                    | -25 to +60<br>-13 to +140               |   |   |

Low-end temperatures require supplemental cooling.

Specifications obtained at sea level using water.

Bath depth includes bezel, height does not include lid.

Thermo Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Thermo Fisher Scientific reserves the right to change specifications without notice.

## General Purpose Water Bath Specifications

|  | GP 02                                  | GP 2S                                  | GP 05                                  | GP 10                                   |
|--|--|--|--|---|
| <b>Reservoir Fluid Control Temperature</b> °C<br>°F          | Ambient to 90<br>Ambient to 194        | Ambient to 100<br>Ambient to 212       |  |   |
| <b>Set Point Temperature Range</b> °C<br>°F                  | +5 to 100<br>+9 to 212                 |  |  |   |
| <b>Ambient Temperature Range</b> °C<br>°F                    | +15 to 45<br>+59 to 113                |  |  |   |
| <b>Stability @ 37°C</b>                                      | ±0.1                                   |  |  |   |
| <b>Uniformity @ 37°C</b>                                     | ±0.2                                   |  |  |   |
| <b>Heating Output</b> watts                                  | 200                                    | 300                                    | 300                                    | 800                                     |
| <b>Bath Volume</b> liters                                    | 2                                      | 2                                      | 5                                      | 10                                      |
| <b>Overall Bath Dimensions</b><br>(L x W x H) cm<br>inches   | 23.0 x 19.9 x 23.3<br>9.1 x 7.8 x 9.2  | 24.6 x 35.5 x 23.2<br>9.7 x 14.0 x 9.1 | 24.6 x 35.5 x 23.2<br>9.7 x 14.0 x 9.1 | 39.3 x 38.3 x 23.3<br>15.5 x 15.1 x 9.2 |
| <b>Bath Work Area Dimensions</b><br>(L x W x H) cm<br>inches | 13.8 x 15.5 x 15.0<br>5.4 x 6.1 x 5.9  | 15.3 x 30.0 x 6.5<br>6.0 x 11.8 x 2.6  | 15.4 x 30.0 x 15.0<br>6.1 x 11.8 x 5.9 | 30.1 x 33.0 x 15.0<br>11.9 x 13.0 x 5.9 |
| <b>Approximate Weight</b> kg<br>lb                           | 3.5<br>7                               | 4<br>9                                 | 4.5<br>10                              | 7.5<br>16                               |
| <b>Electrical Requirements</b> (VAC/Hz)<br>(Voltage ±10%)    | 100-115/50-60<br>or<br>200-230/50-60   |  |  |   |
| <b>Compliance</b>  | CE      RoHS<br>UL      WEEE           |  |  |   |
| <b>Maximum Relative Humidity</b><br>(Non Condensing)         | 80% (up to 31°C)<br>80% (up to 88°F)   |  |  |   |
| <b>Operating Altitude</b> meters<br>feet                     | Sea Level to 2000<br>Sea Level to 6560 |  |  |   |
| <b>Overvoltage Category</b>                                  | II                                     |  |  |   |
| <b>Pollution Degree</b>                                      | 2                                      |  |  |   |
| <b>Storage Temperature Range</b> °C<br>°F                    | -25 to +60<br>-13 to +140              |  |  |   |

Low-end temperatures require supplemental cooling.

Specifications obtained at sea level using water.

Thermal beads may be used instead of water but they will degrade the bath's uniformity and stability.

Bath depth includes bezel, height does not include lid.

Thermo Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Thermo Fisher Scientific reserves the right to change specifications without notice.

## General Purpose Water Bath Specifications

|  | GP 20                                   | GP 28                                    | GP 15D                                  |
|--|---|--|---|
| <b>Reservoir Fluid Control Temperature</b> °C<br>°F          | Ambient to 100<br>Ambient to 212        |  |   |
| <b>Set Point Temperature Range</b> °C<br>°F                  | +5 to 100<br>+9 to 212                  |  |   |
| <b>Ambient Temperature Range</b> °C<br>°F                    | +15 to 45<br>+59 to 113                 |  |   |
| <b>Stability @ 37°C</b>                                      | ±0.1                                    |  |   |
| <b>Uniformity @ 37°C</b>                                     | ±0.2                                    |  |   |
| <b>Heating Output</b> watts                                  | 1200                                    | 1200                                     | 300 and 800                             |
| <b>Bath Volume</b> liters                                    | 20                                      | 28                                       | 5 and 10                                |
| <b>Overall Bath Dimensions</b><br>(L x W x H) cm<br>inches   | 39.2 x 55.5 x 23.3<br>15.4 x 21.8 x 9.2 | 39.2 x 55.5 x 28.2<br>15.4 x 21.8 x 11.1 | 39.2 x 58.7 x 23.3<br>15.4 x 23.1 x 9.2 |
| <b>Bath Work Area Dimensions</b><br>(L x W x H) cm<br>inches | 29.7 x 50.0 x 15.0<br>11.7 x 19.7 x 5.9 | 29.7 x 50.0 x 20.0<br>11.7 x 19.7 x 7.9  | See GP 05 and GP 10                     |
| <b>Approximate Weight</b> kg<br>lb                           | 10<br>22                                | 12<br>26                                 | TBD<br>TBD                              |
| <b>Electrical Requirements</b> (VAC/Hz)<br>(Voltage ±10%)    | 100-115/50-60<br>or<br>200-230/50-60    |  |   |
| <b>Compliance</b>  | CE RoHS<br>UL WEEE                      |  |   |
| <b>Maximum Relative Humidity</b><br>(Non Condensing)         | 80% (up to 31°C)<br>80% (up to 88°F)    |  |   |
| <b>Operating Altitude</b> meters<br>feet                     | Sea Level to 2000<br>Sea Level to 6560  |  |   |
| <b>Overvoltage Category</b>                                  | II                                      |  |   |
| <b>Pollution Degree</b>                                      | 2                                       |  |   |
| <b>Storage Temperature Range</b> °C<br>°F                    | -25 to +60<br>-13 to +140               |  |   |

Low-end temperatures require supplemental cooling.

Specifications obtained at sea level using water.

Thermal beads may be used instead of water but they will degrade the bath's uniformity and stability.

Bath depth includes bezel, height does not include lid.

Thermo Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Thermo Fisher Scientific reserves the right to change specifications without notice.

## Coliform Bath Specifications

|  | COL 19  | COL 35                                  |
|--|---|---|
| <b>Reservoir Fluid Set Point Temperatures</b> °C<br>°F       | 35.0, 41.5, 44.5, 45.5<br>95.0, 106.7, 112.1, 113.9 |   |
| <b>Ambient Temperature Range</b> °C<br>°F                    | +15 to 45<br>+59 to 113                             |   |
| <b>Stability @ 37°C</b>                                      | ±0.1  |   |
| <b>Uniformity @ 37°C</b>                                     | ±0.05   |   |
| <b>Heating Output</b> watts                                  | 1200  | 1500                                    |
| <b>Bath Volume</b> liters                                    | 19  | 35                                      |
| <b>Overall Bath Dimensions</b><br>(L x W x H) cm<br>inches   | 39.4 x 63.2 x 24.9<br>15.5 x 24.9 x 9.8             | 39.4 x 93.8 x 24.9<br>15.5 x 36.9 x 9.8 |
| <b>Bath Work Area Dimensions</b><br>(L x W x H) cm<br>inches | 30.5 x 38.7 x 19.2<br>12.0 x 15.3 x 7.6             | 30.5 x 69.2 x 19.2<br>12.0 x 27.3 x 7.6 |
| <b>Approximate Weight</b> kg<br>lb                           | 16<br>35  | 27<br>60                                |
| <b>Electrical Requirements</b> (VAC/Hz)<br>(Voltage ±10%)    | 100-115/50-60<br>or<br>200-230/50-60                |   |
| <b>Pump Pressure</b> mbar                                    | 262   |   |
| <b>Pump Flow</b> lpm   | 10.6  |   |
| <b>Compliance</b>  | CE      RoHS<br>UL      WEEE                        |   |
| <b>Maximum Relative Humidity</b><br>(Non Condensing)         | 80% (up to 31°C)<br>80% (up to 88°F)                |   |
| <b>Operating Altitude</b> meters<br>feet                     | Sea Level to 2000<br>Sea Level to 6560              |   |
| <b>Overvoltage Category</b>                                  | II  |   |
| <b>Pollution Degree</b>                                      | 2   |   |
| <b>Storage Temperature Range</b> °C<br>°F                    | -25 to +60<br>-13 to +140                           |   |

Specifications obtained at sea level using water.

Bath depth includes bezel, height does not include lid.

Thermo Fisher Scientific takes no responsibility for damages caused by the selection of an unapproved fluids.

Thermo Fisher Scientific reserves the right to change specifications without notice.

# 3

## Chapter 3 Installation

### Bath Installation

The bath is designed for continuous operation and for indoor use.

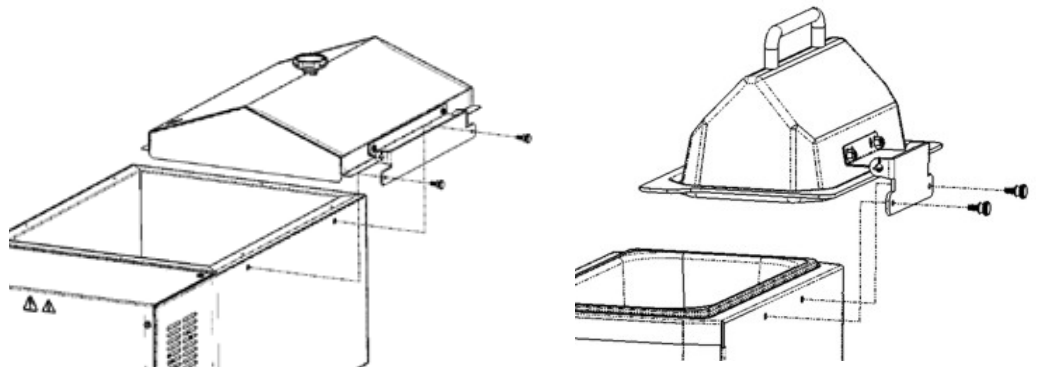
**Never place the bath in a location where excessive heat, moisture or corrosive materials are present.**

**For shaking water baths, the installation platform must be capable of withstanding the shaking generated when the bath is in operation.**

### Ventilation

No special ventilation clearances are required.

### Cover Installation

**Figure 3-1** – Shaking and Circulating Bath Cover Installation, GP Bath Cover Installation

## Electrical Requirements

Refer to the bath's nameplate for specific electrical requirements.

The bath is intended for use on a dedicated outlet.

**The bath construction provides protection against the risk of electrical shock by grounding appropriate metal parts. The protection will not function unless the power cord is connected to a properly grounded outlet. It is the user's responsibility to assure a proper ground connection is provided.**

The circuit protector on the rear of the bath is designed to protect the bath's internal components.

**Note** If the circuit protector activates allow the temperature to cool before resetting. Restart the bath. Contact us if the circuit protector activates again.

**The bath's electrical power cord is the disconnecting device, it must be easily accessible at all times.**

**Ensure the cord does not come in contact with the reservoir contents.**

**Note** Before inserting the electrical cord into the bath connection ensure the circuit protector is in the **O** (off) position.

Once the cord is connected to the bath, connect the other end to the main power source.

## Approved Fluids

The only approved fluids for all Precision baths are:

### Filtered/Single Distilled Water

#### De-ionized Water <sup>1</sup>

<sup>1</sup> For applications requiring resistivity greater than 1 MΩ-cm or maintaining resistivity levels greater than 1 MΩ-cm please call and speak to an applications engineer for additional information.

Thermo Fisher Scientific takes no responsibility for damages caused by using an unapproved fluid.

**When using water above 80°C closely monitor the fluid level, frequent top-offs will be required. It will also create steam.**

### Filtered/Single Distilled Water

Filtered drinking water and single distilled water are good choices for heated water baths. The filtering/distilling process removes microorganisms that are known to cause biological fouling in water baths. The filtering/distilling process also removes minerals and harmful particulates from water. Microorganisms, minerals and particulates can lead to deposits and or scaling that eventually promotes corrosion if not removed or treated.

It is also recommended to regularly drain and replace fluid especially when higher temperatures are utilized (above 40°C). Draining and drying the bath after each use is recommended as leaving standing water for any period of time is known to lead to a buildup of biologic growth.

### Tap Water

Usage of tap water may not cause any adverse affects on the equipment in the short term, but in the long term problems may arise due to precipitation of minerals and or development of algae and bacteria. Please refer to “Process Water Quality and Standard” recommendations in this chapter for guidelines on water usage.

### Chlorine

Draining and replacing “used” water is the best preventative maintenance for the cooling/heating equipment. Biological organisms will develop in water that is not replaced or treated regularly. A mild treatment using an algaecide can help extend the useful life of water.

To help alleviate bacteria or algae growth Thermo Fisher Scientific recommends the use of small amounts of chlorine. The usage of chlorine needs to be monitored over time in order to prevent the formation of corrosion.

The duration of time that chlorine remains in solution depends on factors such as water temperature, pH and availability of direct sunlight. We recommend maintaining chlorine levels

at 1 to 5 ppm free chlorine. PPM levels can be monitored using standard water quality test strips. For best results, maintain the pH of the fluid between 6.5 and 7.5.

Do not add additional chlorine without first determining the concentration ratio that already exists in the fluid supply. Corrosion and degradation of the circulation components can result from concentration ratios that are too high. Contact our customer support for additional information.

### **De-ionized Water**

De-ionized water is water that has had its mineral ions removed using ion exchange resins. The purpose of this process is to remove the ions that allow electrical current to flow more easily through water. This helps to prevent electrical leaks to ground through the recirculating fluid. De-ionized water is in an unbalanced state and will leach the missing ions from the materials it comes in contact with. The aggressive nature of this leaching can cause pitting on metal surfaces. Note that the de-ionizing process does not remove microorganisms. Because of this, we recommend de-ionized water only with applications that have it as a specified requirement. In any case, only de-ionized water with 1 MΩ-cm resistivity, or less, is recommended.

### **Recommended Biocides and Inhibitors**

Thermo Fisher Scientific offers a biocide and inhibitor package Thermo 200 (Nalco) premixed with five gallons of water or as a kit to be added to water.

**Biocides are corrosive and can cause irreversible eye damage and skin burns. They are harmful if inhaled, swallowed or absorbed through the skin. Refer to the manufacturer's most current MSDS.**

### **Thermal Beads (GP Baths only)**

GP baths can use thermal beads. The beads replace water to reduce the risk of contamination.

| Water Quality and Standards |   |                 |
|-----------------------------|---|-----------------|
| Process Fluid               | Permissible (PPM)                           | Desirable (PPM) |
| <b>Microbiologicals</b>     |   |                 |
| (algae, bacteria, fungi)    | 0   | 0               |
| <b>Inorganic Chemicals</b>  |   |                 |
| Calcium                     | <25   | <0.6            |
| Chloride                    | <25   | <100            |
| Copper                      | <13   | <10             |
|                             | 0.020 ppm if fluid in contact with aluminum |                 |
| Iron                        | <0.3  | <0.1            |
| Lead                        | <0.015                                      | 0               |
| Magnesium                   | <12   | <0.1            |
| Manganese                   | <0.05                                       | <0.03           |
| Nitrates/Nitrites           | <10 as N                                    | 0               |
| Potassium                   | <20   | <0.3            |
| Silicate                    | <25   | <1.0            |
| Sodium                      | <20   | <0.3            |
| Sulfate                     | <25   | <1              |
| Hardness                    | <17   | <0.05           |
| Total Dissolved Solids      | <50   | <10             |
| <b>Other Parameters</b>     |   |                 |
| pH                          | 6.5-8.5                                     | 7-8             |
| Resistivity                 | 0.01*                                       | 0.05-0.1*       |

\* MΩ-cm (compensated to 25°C)

Unfavorably high total ionized solids (TIS) can accelerate the rate of corrosion. These contaminants function as electrolytes which increase the potential for galvanic cell corrosion.

Tap water is not normally recommended because the total ionized solids level may be too high. As an example, tap water in the United States averages 171 ppm (of NaCl). The recommended level for use in a water system is between 0.5 to 5.0 ppm (of NaCl).

## Filling Requirements

Ensure the reservoir drain port is closed and that all plumbing connections are securely plumbed or capped. Also ensure any residue is thoroughly removed from the reservoir before filling.

To avoid spilling, place your samples/trays into the bath before filling. Also make allowance for splashing of water from tray oscillation movement.

*Slowly* fill the reservoir. When adding, point the opening of a container away from yourself.

For Shaking Water Baths, fill the reservoir with a minimum of one inch and the maximum level lower than the tray shaft fitting.

For Circulating and Coliform Baths, fill the reservoir with a minimum of two inches, ensuring the suction screen is covered, and a maximum level one inch lower than the upper tank surface.

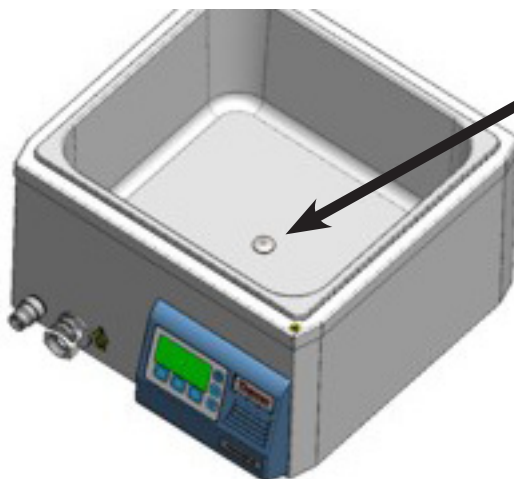
For General Purpose Baths, fill the reservoir with a minimum of one inch and a maximum level one inch lower than the tank upper surface.

**Note** Monitor the fluid level whenever heating the fluid.

To conserve energy, reduce evaporation, and increase temperature control accuracy, use the supplied gable cover. Do not use aluminum foil as a cover, as it may cause corrosion due to an electrochemical reaction.

## Filling with Thermal Beads (GP baths only)

Use the controller to configure the bath before using the thermal bead mode, see Chapter 4 Changing a Setting.

Before filling the reservoir with beads ensure the supplied reservoir plug is securely installed into the drain.

**Note** For safety considerations, enabling the mode automatically derates the bath's heater performance, and it turns off low fluid protection so the beads don't trip the alarms.

This mode also derates the bath's stability and uniformity capability.

**Figure 3-2**– Thermal Beads Reservoir Plug

## Draining

Drain the bath before moving or storing.

Ensure the fluid is at a safe handling temperature, ~40°C or lower. Wear protective clothing and gloves.

- Place a suitable receptacle underneath the drain.
- Slowly turn the drain plug until flow is observed.
- When the flow stops close the drain plug.

### Draining (GP 10, GP 20, GP 28 and GP 15D only)

- Remove the drain hose from its housing and insert the fitting into the drain quick disconnect.
- When the flow stops remove the hose and reinsert it back into the housing.

Figure 3-3– GP Bath Drain Hose

## Accessories

Only use the accessories supplied by Thermo Fisher.

### Drain Kit Accessory

- Before filling the bath remove the drain plug from the bath drain.
- Wrap coupler with Teflon<sup>®</sup> tape, or equivalent, and install into bath drain.
- Insert quick disconnect into hose.
- To start draining snap quick disconnect onto coupler.
- To remove quick disconnect from coupler and stop draining press down on the disconnect's grey release.

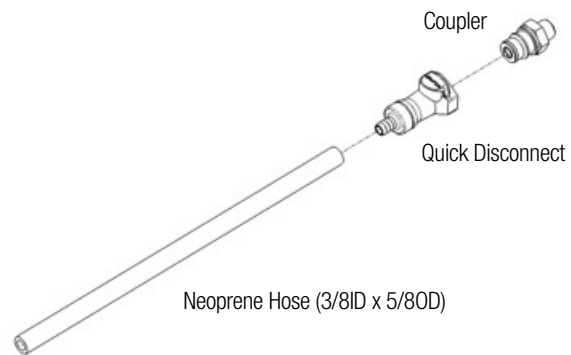
Figure 3-4—Drain Kit Accessory

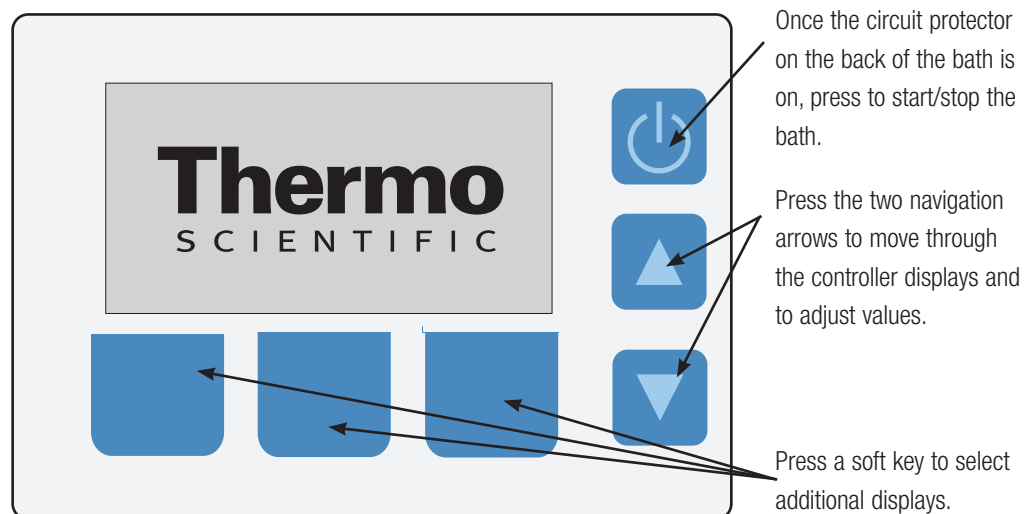
# 4

## Chapter 4 Operation

### Controller

The Thermo Scientific Precision Water Baths have digital controllers that display the bath's reservoir fluid temperature, shaker speed and other bath features.

Indicates the heater is on. The icon flashes when the reservoir fluid temperature is near or at the set point.

Indicates the timer is enabled. The timer is used to turn the bath on, off, or both.

Indicates the audible alarm status, enabled or disabled.

Indicates the shaker mode is enabled, SWB only.

Indicates the pump is running. This icon is always on once the bath is running. Circulating and Coliform Baths only.

## Start Up

Before starting, double check all electrical connections.

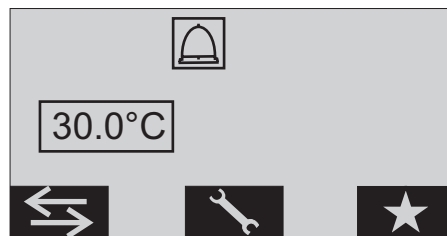
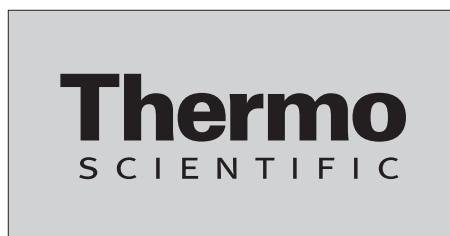
An operating Shaker Water Bath has moving components. Do not insert appendages between the reservoir wall and rack. We recommend stopping the shaker table when loading/removing samples or vessels.

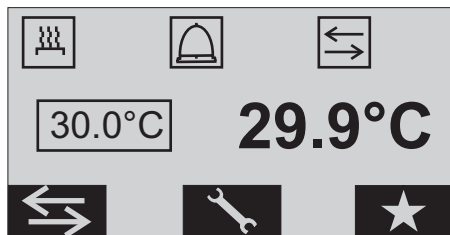
For Shaking Water Baths ensure all glass vessels are secure. Unrestrained vessels may shift and break when shaking is engaged. Flask clips and trays are available from Thermo Fisher.

Do not run the bath until fluid is added to the reservoir. Have extra fluid on hand. If the bath does not start refer to Chapter 6 Troubleshooting.

- Place the circuit protector located on the rear of the bath to the I position.
- The screen will momentarily display **Thermo** and then the home screen appears displaying the set point and the audible alarm status. The set point is the desired reservoir fluid temperature.

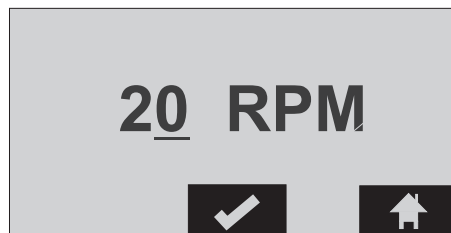
- Press  to start the bath and display the reservoir fluid temperature. If the fluid temperature is below the set point the heater will start.

appears on Shaking Water Baths only.

## Viewing/Changing the Shaker Speed (SWB only)

You can view/change the shaker speed value with the bath running or not. From the home screen press the soft key below  to display:

Press, and hold, the arrow keys to bring up the desired value. The values are 0 and 30 to 200.

Once the desired value is displayed press  or  to save the change and return to the home screen.

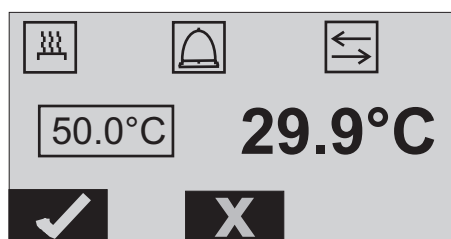
## Changing the Set Point

**Note** This procedure applies to Shaker, Circulating and General Purpose Baths only.

**When operating without a lid, limit the maximum set point to 60°C.**

The set point is the desired reservoir fluid temperature. You can change the set point with the bath running or not. From the home screen press either arrow key to display:

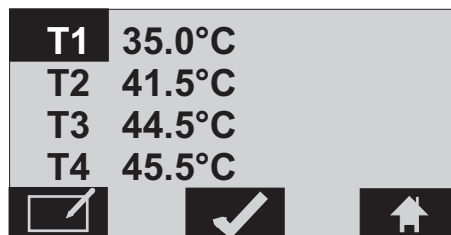
Press, and hold, the arrow keys to bring up the desired set point value. The range is 5°C to 100°C.

Once the desired value is displayed, in this case 50.0°C, press the soft key below  to save the change or press the soft key below  to not save the change. The display returns to the

home screen. **Note** After a 30 second delay if neither key is pressed the display returns to the home screen, any change is *not* saved.

## Selecting a Preset Set Point

The set point is the desired reservoir fluid temperature. You can select the set point with the bath running or not. Or, from the home screen press .

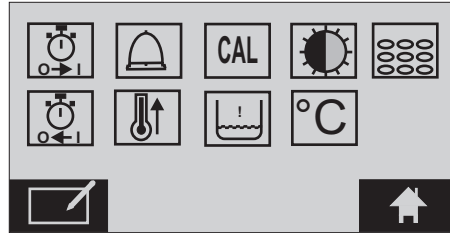
Use the arrow keys to select the desired set point value. Once the desired set point is highlighted press the soft key below  to save the change and return to the home screen.

**Note** After a 30 second delay if neither key is pressed the display returns to the home screen, any change is *not* saved.

**Note** From the home screen you can cycle through the set points by pressing the arrow keys (Coliform baths only).

## Changing a Setting

Press the soft key below  to display the Setting screen.

 appears on General Purpose Baths only.

Use the arrow keys to highlight the desired setting. Press

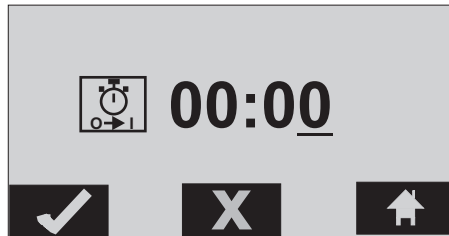
 to make changes.

These icons represent the timer. The top one is used to set the amount of time required, in hours and minutes, before the bath will start. The bottom is used to set how long the bath will run, in hours and minutes, before shutting down. When the bath starts/stops the alarm, if enabled, will beep.

Once the desired timer mode is selected press, and hold, the arrow keys to change the time.

Once the desired value is displayed press  to save the change or  to reset the time to zero. Press  to return to the home screen.

The timer can be used to turn the bath on and then, at a later time, turn the bath off. (Or turn a running

bath off and then back on.) The bath off time must take into account the start time value. For example, if you want the bath to start in 60 minutes and then turn it off 30 minutes later, set the start time to 1:00 and set the off time to 1:30.

**Note** You cannot set a start time if the bath is running or set an off time if the bath is not running.

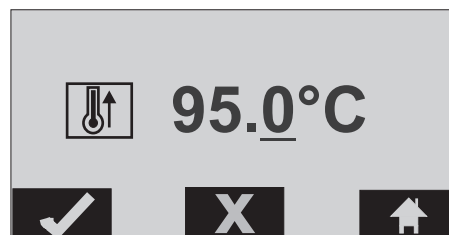
These icons represent the audible alarm status. The top indicates it is enabled, the bottom indicates disabled. Press  to toggle between enabled/disabled. Use the arrow keys to highlight another setting or press

to return to the home screen. Press, and hold, an arrow key to change the temperature. Once the desired value is displayed press  to save the change or  to not save the change.

This icon represents the high temperature alarm. If the reservoir fluid exceeds this temperature the bath will shut down and, if enabled, the alarm will sound. The controller will also have an error display, .

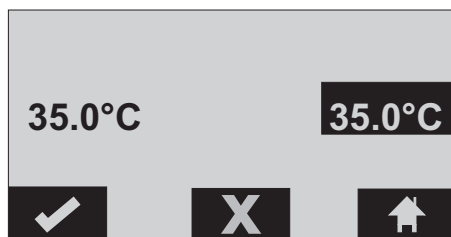
Press  to return to the home screen.

This icon is used to do a reservoir fluid temperature sensor 2-point calibration. The procedure requires a calibrated reference thermometer. Except for Coliform Baths, before the calibration can be performed the fluid temperature must be stabilized at 35°C and then again at 70°C, 35°C and then again at 55°C for Coliform Baths (the sequence does not matter).

**Note** Performing a calibration at any other temperatures will result in an error display .

Once the temperature has stabilized at 35°C, or 70°C (55°C), press  then use the arrow key to change the temperature to match the reference thermometer. Once the desired value is displayed press  to save the change or  to not save the change.

Press  to return to the home screen. Repeat the procedure for the other temperature.

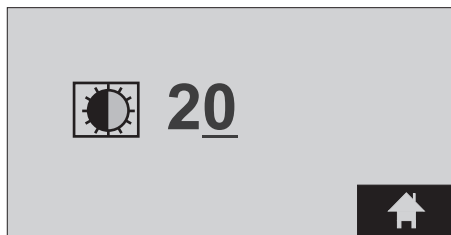
These icons represent the low fluid protection status. The top indicates it is enabled, the bottom is disabled. Press  to toggle between enabled/disabled. Press  to return to the home screen.

When enabled and the reservoir level nears empty the bath will shut down and, if enabled, the alarm will sound. The controller will also have an error display, .

This icon is used to adjust the display's contrast.

Press, and hold, an arrow key to change the contrast. Once the desired contrast is displayed press  to return to the home screen.

These icons are used to select the desired temperature scale, °C or °F. Press  to toggle between °C and °F. Press  to return to the home screen.

GP bath reservoirs can be filled with thermal beads instead of water. These icons are used to select/deselect the reservoir bead mode.

Press  to toggle between on and off. **Note** For safety considerations, the mode automatically derates the bath's heater performance, and it turns off low fluid protection so the beads don't trip the alarms. This mode also derates the bath's stability and uniformity capability.

Press  to return to the home screen.

## Selecting a Preset

Four presets are available. Press the soft key below  to display the Presets screen.

**Note RPM** appears on Shaker Water Baths only.

|   |   |   |
|---|---|---|
| <b>T1</b>   | <b>35.0°C</b>   | <b>30 RPM</b>   |
| <b>T2</b>   | <b>37.0°C</b>   | <b>50 RPM</b>   |
| <b>T3</b>   | <b>50.0°C</b>   | <b>70 RPM</b>   |
| <b>T4</b>   | <b>55.0°C</b>   | <b>90 RPM</b>   |
|  |  |  |

To select the desired preset use the arrow keys to highlight it and then press . The home screen will appear with the desired preset set point value.

## Changing a Preset

**Note** Changing a preset does not automatically select it. Once a change is made you will have to return to the Preset screen and follow the **Selecting a Preset** procedure.

Press the soft key below  to display the Presets screen. **Note RPM** appears on Shaker Water Baths only.

|   |   |   |
|---|---|---|
| <b>T1</b>   | <b>35.0°C</b>   | <b>30 RPM</b>   |
| <b>T2</b>   | <b>37.0°C</b>   | <b>50 RPM</b>   |
| <b>T3</b>   | <b>50.0°C</b>   | <b>70 RPM</b>   |
| <b>T4</b>   | <b>55.0°C</b>   | <b>90 RPM</b>   |
|  |  |  |

To change a preset first highlight it using the arrow keys and then press .

|   |   |   |
|---|---|---|
|  | <b>T1:35.0°C</b>  |   |
| <b>CAL</b>  | <b>RTA1 0.0°C</b>   |   |
|  | <b>30 RPM</b>   |   |
|  |  |  |

Press the arrow keys to highlight the desired preset parameter and press  again. Then press, and hold, the arrow keys to change the value.

Once the desired value is displayed press .

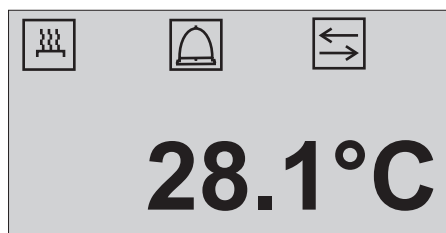
If desired, press the arrow keys to highlight another preset parameter or press  to return to the home screen.

Changing the Real Time Adjustment (CAL RTA) requires a calibrated reference thermometer. If the displayed temperature does not accurately reflect the actual temperature in the reservoir an RTA value is required.

As an example, if the temperature is stabilized and displaying 20°C but a calibrated reference thermometer reads 20.5°C, set the RTA to -0.5°C. After you enter an RTA value allow the display to stabilize before verifying the bath temperature.

**Note** The RTA applies only to the selected preset.

## Temperature Display

With the bath running, and after 30 seconds, if none of the keys are pressed, the controller will display only the reservoir fluid temperature and any active feature. Press any key to return to the previously displayed screen.

## Shut Down

- Press  to stop the bath heating and, if applicable, shaking.

**Note** Pressing and holding  for three seconds stops the heating and powers down the controller display. Press  to turn the display back on.

- For bath shut down place the circuit protector located on the rear of the bath to the **0** position.
- To disconnect the bath remove the electrical plug from its power source.

# 5

## Chapter 5 Preventive Maintenance

**Disconnect the power cord prior to performing any maintenance.**

**Handle the bath with care, sudden jolts or drops can damage its components.**

### **Pump Inlet Filter (Circulating and Coliform baths only)**

The filter on the pump inlet protects the pump from foreign materials. Since applications vary, monitor the filter to establish your cleaning/replacement schedule.

### **Cleaning**

Clean the bath's surface with a soft cloth and warm water only.

Clean the reservoir and built-in components at least every time the bath liquid is changed.

Before cleaning the bath's surfaces, to protect labels, the nameplate, electrical connections, painted and plastic surfaces and to prevent the cleaning agent from entering through any vent openings, mask off all areas except the reservoir.

After time, the circulating bath's stainless steel surfaces may show spots and become tarnished. Normal stainless steel cleaners can be used.

**Do not use scouring powder or any substance containing solvents.**

The inside of the bath must be kept clean in order to ensure a long service life. Quickly remove substances containing acidic or alkaline substances and metal shavings as they could harm the surfaces causing corrosion. If corrosion (e.g., small rust marks) occur in spite of this, cleaning with stainless steel caustic agents has proved to be suitable. Apply these substances according to the manufacturer's recommendations.

## Electrical Power Cord

Ensure any replacement cord is properly rated.

## Optical Disc (SWB only)

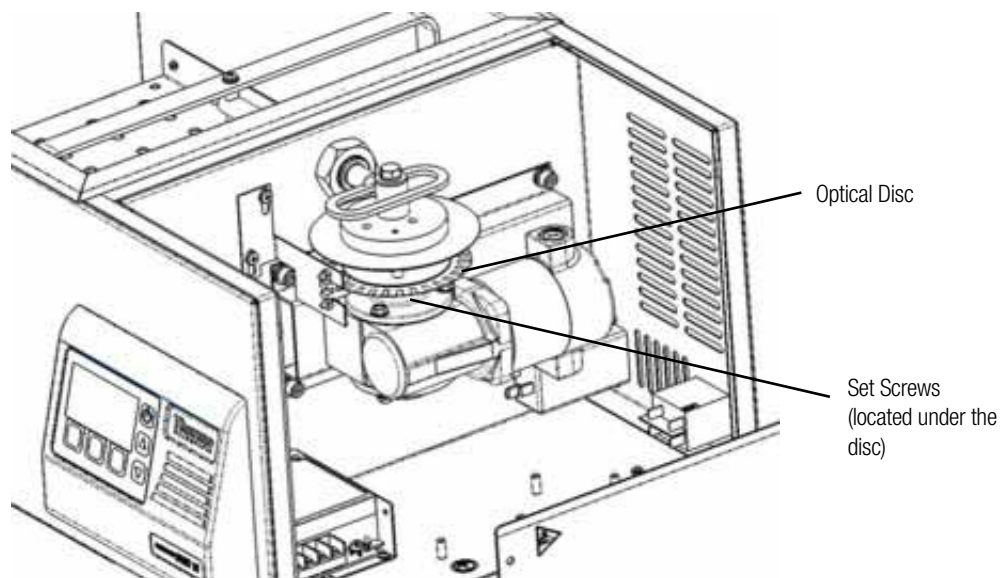
Check the optical disc monthly and verify it is free of foreign material. If not, wipe the disc clean.

1. Remove the four screws, two on top and two on the side securing the top right panel.
2. Swing panel assembly open.

**Do not operate the bath with the panel open.**

3. Rotate disc, by grasping aluminum components of the eccentric only, one turn clockwise and verify that disc does not rub against optical sensor. If required, use the 3/32" Allen to loosen two set screws on the eccentric and reposition the disc.

## Shaker Assembly

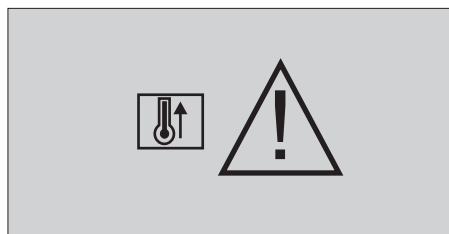
Check the shaker assembly for looseness annually. Contact Fisher Scientific if excessive movement is noted.

# 6

## Chapter 6 Troubleshooting

### Error Displays

An error display indicates an unusual condition. With any message, except **CAL**, the bath will stop heating and, if applicable, the pump will stop. With any message the alarm, if enabled, will sound.

| Display   | Cause/Action   |
|---|--|
|  | <b>HTC Fault</b><br>fixed high temp protection limit exceeded<br>allow bath to cool down<br>restart the bath<br>if the HTC fault cannot be cleared, the bath must be serviced by an authorized Thermo Scientific Temperature Control Service Technician. |
|  | <b>High Temp Fault</b><br>adjustable high temp fault protection limit exceeded<br>check limit setting  |
|  | <b>Low Fluid Level</b><br>low level protection limit exceeded<br>check fluid level<br>check for leaks  |
|  | <b>Open/Shorted Internal Temp Sensor</b><br>the bath must be serviced by an authorized Thermo Scientific Temperature Control Service Technician.   |

| Display   | Cause/Action   |
|---|--|
|  | <p><b>Bad Calibration</b></p> <p>redo calibration</p> <p>ensure calibration set points are 35.0°C and 70.0°C (35.0°C and 55.0°C for Coliform Baths)</p>  |
|  | <p><b>Shaker Malfunction (SWB only)</b></p> <p>ensure there is nothing in the reservoir hampering shaker movement.</p> <p>the bath must be serviced by an authorized Thermo Scientific Temperature Control Service Technician.</p> |

## Checklist

### Bath will not start or shuts down

- Check display for error messages.
- Ensure  wasn't accidentally pressed.
- Ensure the circuit protector is in the on ( I ) position.
- Check the line cord connection to your power supply and at the bath.
- Make sure supply voltage is connected and matches the bath's nameplate rating  $\pm 10\%$ .
- Restart the bath.

### No display

- Pressing  should return the display.
- Cycle the bath's circuit protector.

### Inadequate temperature control

- Verify the set point.
- Low-end temperature set points require supplemental cooling.
- GP 02 baths require insulation, especially in the lid, for operation above 90°C.
- Ensure bath installation complies with the site requirements in Chapter 3.
- Make sure supply voltage matches bath's nameplate rating  $\pm 10\%$ .

Please contact Thermo Fisher Scientific Sales Service and Customer Support if you need any additional information, see inside cover for contact instructions.

# Warranty

Thermo Fisher Scientific warrants products for 24 months from date of shipment according to the following terms.

Any part of the bath manufactured or supplied by Thermo Fisher Scientific and found in the reasonable judgment of Thermo Fisher to be defective in material or workmanship will be repaired at an authorized Thermo Fisher Repair Depot without charge for parts or labor. The bath, including any defective part must be returned to an authorized Thermo Fisher Repair Depot within the warranty period. The expense of returning the bath to the authorized Thermo Fisher Repair Depot for warranty service will be paid for by the buyer. Our responsibility in respect to warranty claims is limited to performing the required repairs or replacements, and no claim of breach of warranty shall be cause for cancellation or rescission of the contract of sales of any bath. With respect to baths that qualify for field service repairs, Thermo Fisher Scientific's responsibility is limited to the component parts necessary for the repair and the labor that is required on site to perform the repair. Any travel labor or mileage charges are the financial responsibility of the buyer.

The buyer shall be responsible for any evaluation or warranty service call (including labor charges) if no defects are found with the Thermo Scientific product.

This warranty does not cover any bath that has been subject to misuse, neglect, or accident. This warranty does not apply to any damage to the bath that is the result of improper installation or maintenance, or to any bath that has been operated or maintained in any way contrary to the operating or maintenance instructions specified in this User Manual. This warranty does not cover any bath that has been altered or modified so as to change its intended use.

In addition, this warranty does not extend to repairs made by the use of parts, accessories, or fluids which are either incompatible with the bath or adversely affect its operation, performance, or durability.

Thermo Fisher Scientific reserves the right to change or improve the design of any bath without assuming any obligation to modify any bath previously manufactured.

THE FOREGOING EXPRESS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO WARRANTIES OR MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

OUR OBLIGATION UNDER THIS WARRANTY IS STRICTLY AND EXCLUSIVELY LIMITED TO THE REPAIR OR REPLACEMENT OF DEFECTIVE COMPONENT PARTS AND Thermo Fisher Scientific DOES NOT ASSUME OR AUTHORIZE ANYONE TO ASSUME FOR IT ANY OTHER OBLIGATION.

Thermo Fisher Scientific ASSUMES NO RESPONSIBILITY FOR INCIDENTAL, CONSEQUENTIAL, OR OTHER DAMAGES INCLUDING, BUT NOT LIMITED TO LOSS OR DAMAGE TO PROPERTY, LOSS OF PROFITS OR REVENUE, LOSS OF THE BATH, LOSS OF TIME, OR INCONVENIENCE.

This warranty applies to chillers sold by Thermo Fisher Scientific. (Refer to the warranty for baths sold by the affiliated marketing company of Thermo Fisher Scientific for any additional terms.) This warranty and all matters arising pursuant to it shall be governed by the law of the State of New Hampshire, United States. All legal actions brought in relation hereto shall be filed in the appropriate state or federal courts in New Hampshire, unless waived by Thermo Fisher Scientific.

[thermoscientific.com](http://thermoscientific.com)

© 2015 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

**Thermo Fisher Scientific**  
81 Wyman Street  
Waltham, MA 02451

**Thermo**  
SCIENTIFIC  
A Thermo Fisher Scientific Brand