

# 4' Protector XL Benchtop Laboratory Hood with 2 services fixtures and 1 electrical duplex, 24" interior depth



View online: <https://www.labconco.com/product/4-protector-xl-benchtop-laboratory-hood-with-2-services-fixtures-and-1-e/3795>

## Overview

---



The Protector XL Benchtop Laboratory Hoods are general chemistry hoods that offer superior design flexibility. They are available in widths from 3 to 16 feet and three depths to meet a variety of installation and application requirements. These hoods incorporate containment-enhancing features allowing them to operate at face velocities as low as 60 fpm and meet the SEFA-1 standard of a low velocity high performance hood.

[\*\*SEE OUR SUSTAINABILITY INITIATIVES\*\*](#)

**Catalog Number:** 111400002

## Specifications

---

- **Weight:** 375.0 lbs
- **Weight metric:** 170.0 kg
- **Dimensions:** 48.0" w x 31.7" d x 59.0" h
- **Dimensions metric:** 121.9 x 80.5 x 149.9 cm
- **Electrical:** 100-115 volts, 50/60Hz, 10 amps
- **Product Subcategory:** General Chemistry
- **Nominal Width:** 4'
- **Blower Requirements:** Remote blower required
- **Conformance:** ANSI Z9.5, ASHRAE 110, ASTM E84, CAN/CSA C22.2, NFPA 45, SEFA 1, UL
- **Electrical Duplexes:** 1
- **Enclosure Height:** 59.0"
- **Exterior Depth:** 31.7" (80.5 cm)
- **Lighting:** LED
- **Service Fixtures:** 2
- **Style:** Benchtop

# Description

---

## Features

- By-pass airflow design
- Eco-Foil Air Foil with aerodynamic Clean-Sweep™ airflow openings
- Cord-Keeper™ Slots on left and right side of air foil
- Glacier white powder-coated steel exterior
- Chemical-resistant, fiberglass-reinforced, composite panel liner and baffle
- Opti-Zone™ Baffle with flame spread index less than 25 per ASTM E84. Baffle is removable for cleaning
- Tempered safety glass vertical-rising sash(es) with cable and pulley(s) and powder-coated aluminum sash handle(s). 8' models have optional split dual sashes. 10', 12' and 16' models have split dual sashes (four on 16' models)
- 37.5" (95.3 cm) high sightline from the work surface and header panel
- Removable front and side panels, and front and interior service access panels for access to plumbing and electrical wiring
- Pre-wired LED lighting with vapor-proof design and ADA-compliant light and blower switches. Lighting rated to 50,000 hours
- Powder-coated stainless steel, 12.8" (32.5 cm) ID exhaust connection(s)
- 10', 12' and 16' models also feature a by-pass block

## Required Accessories (not included)


- Remote blower
- Ductwork
- Work Surface
- Base cabinet or stand

## Optional Accessories (not included)


- Service Fixture Kits (some models include service fixtures)
- Electrical Duplex Kits (some models include electrical duplex)
- Guardian Airflow Monitor Kits
- Distillation Grid Kits
- Snuffer™ Fire Extinguisher
- Ceiling Enclosure Kits
- Rear Finish Panel Kits

# For Product Assistance

---

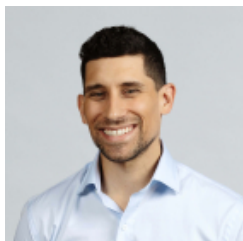

[Ryan Caton](#)  
Sales Representative

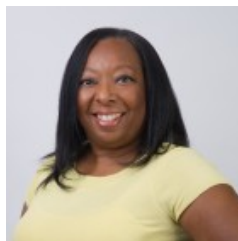

[Eleanor Watkins](#)  
Customer Service Representative



Labconco Corporation  
8811 Prospect Avenue  
Kansas City, MO 64132-2696  
Phone: (816) 333-8811

