



IsoPure™ Consumables


Accuris offers magnetic tip combs, deep well plates, and 6-well sample strips designed to fit the IsoPure Mini. Produced using virgin polypropylene and packaged in a sterile environment, these consumables are guaranteed to be contaminant-free. The 8-place magnetic tip combs lock into position inside the instruments processing chamber. There are 2 options for sample vessels: 96 deep-well plates can be used for processing up to 16 samples. For 8 or fewer samples, 6-well sample strips can be used, and these are loaded onto the included aluminum rack. These sample strips reduce waste and lower consumable costs when only a few samples need to be processed.

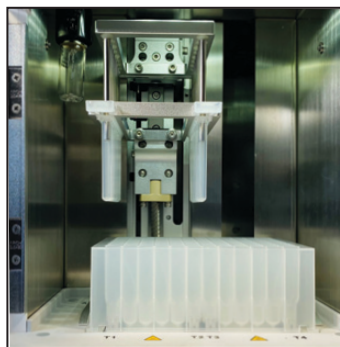

Deep-well plate for processing up to 16 samples

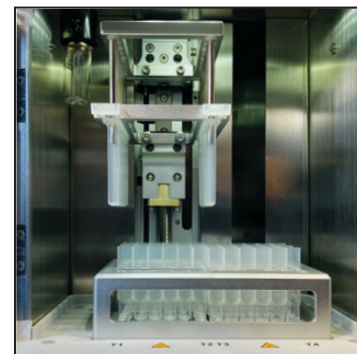

6-well sample strips loaded in aluminum rack

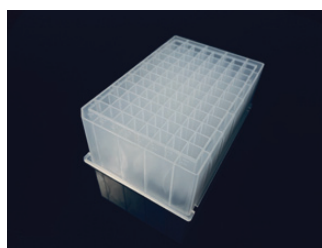


**AP1096-DWP
96 deep well plate**



**AP1016-ASR
Aluminum rack**



**AP1016-6W
6-well sample strip**

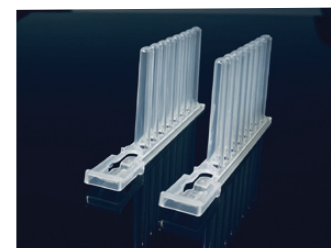

**AP1016-MC
8-place tip comb**

Specifications

| | |
|-----------------------------|--|
| Display | 4.3 inch-color touch screen |
| Capacity / sample volume | 16 samples / 50µL to 1000µL |
| Magnetic bead recovery rate | >95% |
| Lysis temperature range | Ambient +5°C to 120°C |
| Elution temperature range | Ambient +5°C to 120°C |
| Mixing speed | Variable, 10 levels |
| Programming control | Bluetooth-enabled android device(not included) |
| Program memory | 100 programs (internal memory) |
| Communication ports | USB for mouse, file storage & transfer |
| Internal lighting | LED |
| Contamination control | UV lamp (3W, 253.7nm UVC wavelength) |
| Program Transfer | Bluetooth, QR scanner, USB flash drive |
| Dimensions | 7.9 x 10.2 x 11.8 in. (20 x 26 x 30 cm) |
| Weight | 7.0 kg / 15.4 lbs |
| Electrical input | Universal 120 to 240VAC |
| Warranty | 2 years |

Ordering Information

| Cat No. | Description |
|------------|---|
| AP1016* | IsoPure™ Mini, Automated Purification System, 120V - 240V input |
| AP1096-DWP | 96 deep well plates for IsoPure™, pack of 50 |
| AP1016-MC | 8-place magnetic tip comb, pack of 50 |
| AP1016-6W | 6-well sample strips, pack of 100 |
| AP1016-TAB | Android tablet preloaded with IsoPure™ app |

* Use part number AP1016-E to indicate 240V input and specify required plug type.

Also available from Accuris Instruments



IsoPure™ 96

- Process up to 96 samples in deep-well plates
- Compatible with all popular magnetic bead kits
- Touch screen panel for on-board programming & operation
- 8 position carousel with 2 heating blocks
- Sample volume 50µL-1000µL