


# Accuris Balance Accessories

## Accuris Anti Static Ionizer

Static charge can accumulate on many types of samples as well as weighing containers, affecting mass determination on an analytical balance. The Accuris ASI-100 produces and emits negatively and positively charged ions that will neutralize static electricity on samples and containers. A blower can be turned on for rapid static discharge, or turned off for work with powders or light-weight materials.

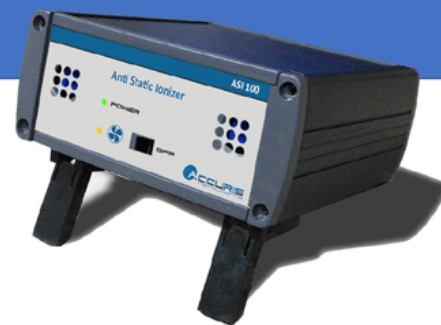

Compact Ionizer  
Order No. W3120

## Portable Thermal Printer

The W3130 Thermal Printer offers a small footprint, modern design, and is compatible with Accuris series *Dx* and *Tx* Analytical Balances. Connection is easy with the included cables, and weight values and other critical info can be printed quickly and clearly. The W3130 printer connects to AC power (100V to 240V) with an included power adapter, and also features an included rechargeable Li-Ion battery.



Thermal Printer  
Order No. W3130



## Calibration Weights

Laboratory balances, from economy to expensive analytical models, require routine calibration and linearity confirmation. Accuris™ offers a selection of individual weights and sets compatible with different types of balances.

Our Premium Class E2 Weight Set is ideal for use with analytical balances. Precision manufactured from high purity, 316 stainless steel ( $7.96\text{kg/dm}^3$ ), these weights are corrosion resistant with magnetic susceptibility  $<0.005$ . The set includes 6 weights, 1x100mg, 1x10g 2x20g, 1x100g, for checking repeatability, corner load, linearity and span. The weights are supplied in a protective case with tweezers, a cleaning cloth and a certificate of verification.

### Ordering Information

W3101A-120	Analytical Balance, series Dx, internal calibration, Graphical Display, 120gx0.0001g
W3101A-220	Analytical Balance, series Dx, internal calibration, Graphical Display, 220gx0.0001g
W3002A-120	Analytical Balance, series Tx, internal calibration, Touch Screen, 62gx0.00001g /120gx0.0001g
W3102A-220	Analytical Balance, series Tx, internal calibration, Touch Screen, 220gx0.0001g
W3120	Compact Anti Static Ionizer
W3130	Thermal Printer for Accuris Series Dx and Tx Balances
W1101-6-100	Calibration weight set, Class E2, for Analytical Balances, 6 pcs: 100mg to 100g



A division of Benchmark Scientific  
PO Box 709 Edison, NJ  
Phone: 908-769-5555  
[www.accuris-usa.com](http://www.accuris-usa.com)