

150 ml Complete Fast-Freeze Flask



View online: <https://www.labconco.com/product/150-ml-complete-fast-freeze-flask/3191>



Overview

Fast-Freeze Flasks attach with adapters to valve ports on Labconco Drying Chambers and Manifolds. They are available in clear or amber-colored glass. A complete Fast-Freeze Flask consists of a rubber top, glass bottom and a supply of filter paper.

Tops, bottoms and filter paper are available separately as replacement components. Adapters are available in borosilicate glass or stainless steel. Fast-Freeze Flasks and Adapters are compatible with any brand freeze dryer that has 1/2" or 3/4" diameter valve ports.

Glass flasks withstand maximum temperature of 490° C. Flask toppers withstand a maximum temperature of 450° C.

Complete Flask consists of 7542400 Flask Bottom, 7544200 Flask Top and a supply of filter paper.

Catalog Number: 7540400


Specifications

- **Weight:** 1.0 lbs
- **Weight metric:** 1.0 kg
- **Dimensions metric:** Flask Bottom: 85mm high, 59.2mm dia.
- **Flask Top Adapter Diameter:** 3/4"
- **Volume:** 150 ml

For Product Assistance



[KeiShaun Johnson](#)
Customer Service Representative

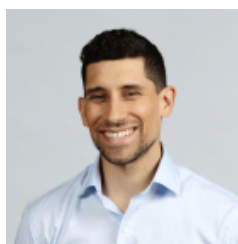

[Ryan Caton](#)
Sales Representative



Labconco Corporation
8811 Prospect Avenue
Kansas City, MO 64132-2696
Phone: (816) 333-8811

