



powder coated, heavy-duty steel body is easy to clean. Temperature is displayed on an easy-to-read LED readout for Models 26025 and 26025-1.


MODEL 26005

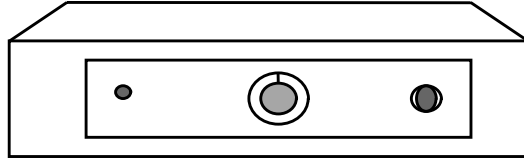

MODEL 26020

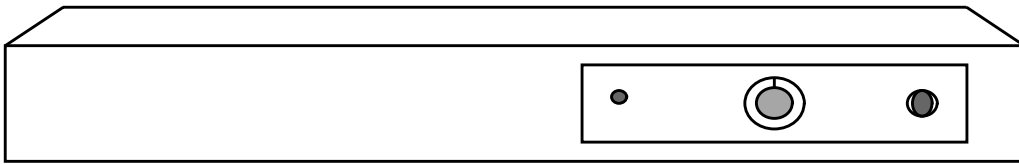

MODEL 26025

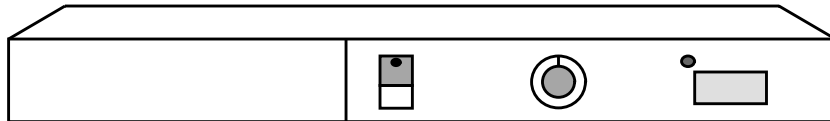

SECTION 3

## SPECIFICATIONS

POWER REQUIREMENTS:	26005:	120 VAC,	50/60 Hz, 160 Watts
	26007:	240 VAC,	50/60 Hz, 160 Watts
	26020:	120 VAC,	50/60 Hz, 160 Watts
	26022:	240 VAC,	50/60 Hz, 160 Watts
	26025:	120 VAC,	50/60 Hz, 160 Watts
	26025-1:	240 VAC,	50/60 Hz, 160 Watts

TEMPERATURE RANGE: Slightly above ambient to 65°C.

SLIDE CAPACITY:	MODELS	CAPACITY	SLIDE SIZE
	26005, 26007	65 32	1" x 3" or 2" x 3"
	26020, 26022	66 33	1" x 3" or 2" x 3"
	26025, 26025-1	66 33	1" x 3" or 2" x 3"

DIMENSIONS:	MODELS	DIMENSIONS
	26005, 26007	14" L x 14" W x 3½" H (36 x 36 x 9 cm)
	26020, 26022	25" L x 8" W x 3½" H (64 x 20 x 9 cm)
	26025, 26025	25" L x 8" W x 4½" H (64 x 20 x 11 cm)

NOTE: CE MODELS HAVE SAME SPECS. AS CORRESPONDING 240 VAC UNITS.

*UNIT'S ENVIRONMENTAL OPERATING CONDITIONS:*

POLLUTION DEGREE:	2
INSTALLATION CATEGORY:	II
ALTITUDE:	2000 Meters MSL (Mean Sea Level)
HUMIDITY:	80% maximum, non-condensing
ELECTRICAL SUPPLY:	120VAC or 240VAC
VOLTAGE TOLERANCE:	±10% of normal rated line
TEMPERATURE:	15°C to 40°C
PRODUCT USAGE:	This product is intended for use indoors <b>only</b>

SECTION 4

## INSTALLATION

### √SHIPPING CARTON:

This should be inspected upon delivery. When received, carefully examine for any shipping damage before unpacking. If damage is discovered, the delivering carrier should both specify and sign for the damage on your copy of the delivery receipt.

Open the carton carefully making certain that all parts are accounted for before packaging materials are discarded—after unpacking, if damage is found, promptly report it to the carrier and request a damage inspection promptly.