

VWR® ANALOG ROCKING PLATFORM SHAKER



VWR® Analog Rocking Platform Shaker

- Variable control for speed, tilt and time
- 16 lb capacity
- Two-tier model doubles workable capacity
- 5 year warranty on parts and labor

VWR® Analog Rocking Shakers are an easy and economical option for all of your rocking needs. Ideal for cell culture and blotting applications and is designed to be used in a variety of environmental conditions.

Operating Features:

Low profile design: Takes up less bench space and fits in most hoods and incubators. Two-tier option increases the capacity with the same footprint and provides a 3.5" (8.9cm) clearance between platforms. Cast aluminum base offers durability and added stability.

Microprocessor control: Provides tilt adjustment which allows user to easily change rocking angle from 0 to 15° while unit is operating. Smooth speed control with low speed rocking motion.

Independent control knobs for speed, tilt, and time, allow for easy adjustments.

Safety Features:

Overload protection: Audible signal will activate when system detects an obstruction or overload of the tray.

Speed ramping feature: Slowly increases speed to desired set-point to avoid splashing.

Timer will automatically stop rocking motion when timer reaches zero.

Spill-resistant design: Channels fluids away from internal components.

Operating Conditions:


Unit can be run in cold rooms, incubators, and CO₂ environments from -10 to 60°C (14 to 140°F), maximum 80% relative humidity, non-condensing.

Applications:

Staining and destaining gels, hybridization procedures, hematology, and blotting techniques.

Ordering Information:

Units include a 92" (234cm) detachable, 3-wire cord and plug (230V units are supplied with a Euro type plug) Units are also supplied with a 14 x 11" (35.6 x 27.9cm) non-skid rubber mat. Two-tier models include a 2nd tray, non-skid mat and hardware. Optional Dimpled Mat for securing tubes and Stacking Tray for single tier rocker. 5 year limited warranty on parts and labor. TUV listed. 230V units are CE marked.

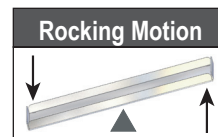

Specifications	
Speed Range	1 to 75rpm*
Tilt Angle	0 to 15° *
Timer	1 minute to 120 minutes
Maximum Weight Capacity	16lbs (7.3kg)**
Tray Material	Aluminum
Tray Dimensions (L x W)	14 x 11" (35.6 x 27.9cm)
Overall Dimensions (L x W x H)	1-Tier 17 x 11 x 5" (43.2 x 27.9 x 12.7cm) 2-Tier 17 x 11 x 9" (43.2 x 27.9 x 22.9cm)
Ship Weight	15.5lbs (7kg)

* Maximum speed/tilt angle may vary with heavy or unbalanced loads.

** Centered on tray.

VWR® Rocking Platform Shaker

Description	Electrical (50/60 Hz)	Cat. No.
One Tier-Rocking Platform Shaker	120V 0.125 amps 15 watts	10127-872
One Tier Rocking Platform Shaker	230V 0.065 amps 15 watts	10127-874
Two-Tier Rocking Platform Shaker	120V 0.125 amps 15 watts	10127-876
Two-Tier Rocking Platform Shaker	230V 0.065 amps 15 watts	10127-724
Dimpled Mat		10127-728
Stacking Tray		10127-726



1.800.932.5000 | vwr.com

Prices and product details are current when published; subject to change without notice. | Certain products may be limited by federal, state, provincial, or local regulations. | VWR makes no claims or warranties concerning sustainable/green products. Any claims concerning sustainable/green products are the sole claims of the manufacturer and not those of VWR International, LLC. All prices are in US dollars unless otherwise noted. Offers valid in US and Canada, void where prohibited by law or company policy, while supplies last. | VWR, the VWR logo and variations on the foregoing are registered (®) or unregistered trademarks and service marks, of VWR International, LLC and its related companies. All other marks referenced are registered by their respective owner(s). | Visit vwr.com to view our privacy policy, trademark owners and additional disclaimers. ©2014 VWR International, LLC. All rights reserved.